

## Bite-Chek® Instructions for Use

**CAUTION:**

Rx Only. These instructions, in whole or in part, are not a substitute for formal training. Appropriate professional education is STRONGLY RECOMMENDED prior to using this device clinically.


**DESCRIPTION:**

Bite-Chek is made up of layers of laminate, ink, paper, medical grade adhesive, and coated film. It is to be used by qualified dental practitioners in dental clinics, hospital, lab, or schools for a dental application. Articulation film is a non-invasive device that can be used on children, adults, adolescents, or animals. Product claims to mark occlusal surfaces or contact points. It is offered as not sterilized; however, it is packaged in a controlled environment. Product is single-use and should NOT be re-processed and/or re-used.

**INDICATIONS**

Bite-Chek is an articulation film product made to mark the occlusal surface of the mouth. Marking surfaces can include zirconia, porcelain, composite, gold, eMax, or stainless steel.

**CONTRAINDICATIONS TO USE**

Use of Bite-Chek is contraindicated on any patient who is allergic to any of the components of the product. Do not reuse. The Bite-Chek are single-use. 

**CLINICAL PRECAUTIONS AND WARNINGS:**

- a) Carefully read package labels to ensure use of the appropriate device.
- b) Failure to follow instructions may cause procedural delays or patient or user injury.
- c) Prior to use, carefully inspect the product for signs of damage and/or deterioration.
- d) Discard any damaged Bite-Chek immediately.
- e) Bite-Chek are for SINGLE-PATIENT-USE ONLY, in a dental setting.
- f) Discard immediately after use.
- g) Always keep track of Lot Numbers of Bite-Chek to ensure traceability.
- h) Used by a trained dental practitioner, assistant, or hygienist.


**CLINICAL USE:****Step #1**

- Dry the tooth surface.



**Step #2**

- Position so film is on the teeth facing the center of the mouth.

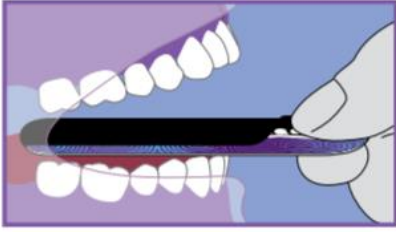

**Step #3**

- Check occlusion



**STORAGE**

Bite-Chek should be stored in a dry, dust-free environment, outside of direct sunlight. Improper storage conditions may cause malfunction of the product.

**TRACEABILITY**

Each package includes **Lot number** LOT on its label.

This number must be quoted in any correspondence regarding the product.

**SYMBOLS:**

	Manufacturer	Indicates the medical device manufacturer.		Do not reuse	Indicates a medical device that is intended for one use, or for use on a single patient during a single procedure.
<span style="border: 1px solid black; padding: 2px;">LOT</span>	Lot Code	Indicates the manufacturer's batch code so that the batch or lot can be identified.		Consult instructions for use	Indicates the need for the user to consult the instructions for use.
	CE marking	Signifies European technical conformity.		Do not use if package is open or damaged	Indicates a medical device that should not be used if the package has been damaged or opened.
<span style="border: 1px solid black; padding: 2px;">EC REP</span>	Authorized European representative	Indicates the Authorized representative in the European Community.			

CONTACT INFORMATION:



Microcopy Dental  
3120 Moon Station Rd. NW  
Kennesaw, GA 30144



DENTEQ Medical Technologies  
Hafenstrasse 12 76344,  
Egenstien-Leopoldshafen, Germany  
Fax: +49-7247 944843



REVISION HISTORY:

MCD-IFU-005 Rev: 1  
Date of issue: 13Nov2018

**OBSOLETE**